

Inochi Care

Healing The
Non-Healing Wounds

InoHeal-Mini

PRECISE
PORTABLE
PROVEN

InoHeal-Mini

Single use, ultra portable Negative Pressure Wound Therapy.

Inoheal mini is a single use Negative Pressure Wound Therapy (sNPWT) system. It is a portable, canister-free medical device that applies gentle suction to a wound to remove fluid, reduce infection risk, and promotes wound healing.



-40 | - 80 | -120 mmHg



Up to 28-day runtime





inochi care
innovating healthcare

About Us

Inochi Care is a healthtech start up developing innovative solutions for wound care and management.

Inochi Care is an IP and tech driven company with focus on wound care products, products that are designed and engineered for scale and adoption across healthcare sectors and different markets.

Meet InoHeal-Mini from Inochi Care



Portable



Canister-free



Designed for
convenient wound
management



InoHeal-Mini is a single-use Negative Pressure Wound Therapy (sNPWT) system that applies controlled negative pressure to the wound.

It is designed to help:

- ✓ Remove excess wound exudate and fluid
- ✓ Maintain controlled negative pressure at the wound site
- ✓ Support the wound healing process
- ✓ Provide convenient, portable negative pressure therapy

InoHeal-Mini



3 Pressure Levels

Choose -40, -80 or -120 mmHg according to the wound and prescribed therapy.



Canister-Free Design

A compact system without the need of a separate fluid collection canister.



Replaceable AA Batteries

Replace batteries according to usage without replacing the entire device for up to four weeks.



Feedback Alarms

Alerts for leaks, blockages, low battery and therapy status.



Dressings

Matters.

InoHeal dressings are available in different sizes and feature an active functional layer designed to support the wound healing process.



Designed for different wound requirements. The dressing connects directly to InoHeal-Mini, creating a simple, compact NPWT solution.

Available Size	Dimension
Small	15 x 10 cm
Medium	15 x 30 cm
Large	10X30 cm

Key Features

- ✓ Adjustable Pressure
- ✓ Up to 30-Day Runtime
- ✓ Intelligent Feedback
- ✓ Replaceable Batteries
- ✓ Portable Design
- ✓ Cost-Efficient

A cost-efficient alternative to imported NPWT systems.

Key Benefits of Single-use Negative Pressure Wound Therapy

- ✓ Improves healing rates
- ✓ Reduces surgical site infections
- ✓ Shortens median healing times compared to standard dressings




How to Use

Follow the five simple steps to start the therapy.


1. Apply the Dressing

Apply dressing over the clean wound and ensure there are no folds.

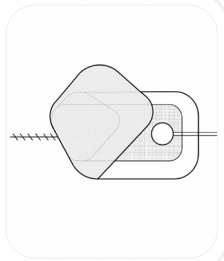


2. Insert batteries

Insert two (2) AA batteries into the battery slot.

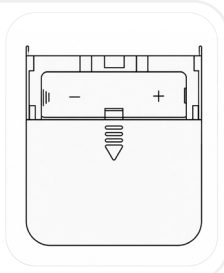


3. Connect tube

Connect the dressing tube to the device.

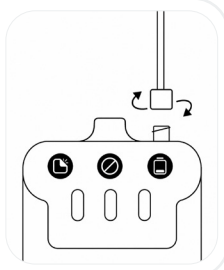


4. Set pressure

Set the pressure between -40, -80, and -120 mmHg with the pressure button.

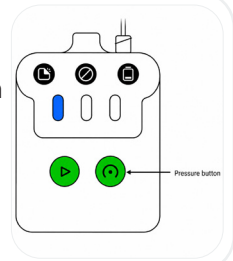

5. Start the therapy

Press the start button to start the therapy

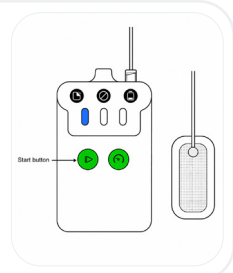

InoHeal-Mini

at a Glance

Specification	Details
Therapy	Single use Negative Pressure Wound Therapy
Device Type	Portable, canister-free
Pressure Levels	40/80/120 mmHg
Battery	2 x AA batteries
Runtime	Up to 30 days
Device Dimensions	10 × 6.5 × 3.5 cm
Dressing Options	Small / Medium / Large
Feedback	Therapy / Leak / Block / Low Battery

Package Contents



- ✓ 1 × InoHeal-Mini device
- ✓ 3 × InoHeal Dressings
- ✓ 2 × AA batteries
- ✓ 1 × User Manual

A compact approach to negative pressure wound therapy.



inochi care
innovating healthcare

For more information, contact

 +91 11404 15913  contact@inochihealthcare.com

 www.inochihealthcare.com